

Alert-it Alarm Message Responses



Name of user **Alarm:**

URGENT = Life critical alarm: Attend immediately

HELP = Comfort critical alarm: Attend with speed

ASSIST = Non critical alarm

FAULT = Equipment needs servicing


Special Fault Alarms

RF FAIL: Correct immediately as radio contact has been lost with the user's equipment

FAULT30: The user's battery needs charging or replacing, when convenient

FAULT31: The user's battery is almost exhausted and needs charging re replacing immediately

Action:

1. The pager sound can be muted by pressing 
2. Attend the user.
3. Reset the alarm on the user's equipment (see examples below)



Press to reset



Press for 1 second until LED stops flashing (and beep heard on some units)

Press *Disarm* to cancel



Full adjustment details are found in handbooks available on:
www.alert-it.co.uk/handbooks/
Or by phoning Alert-it

Basic User Pager Information

Typical screens in normal operation




or



Press 3 sec to enter setup



Press to silence alarms

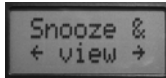


Scroll to menu item required eg.



Select item

(shows a list of all alarm units)

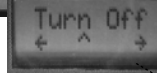

Turn off may require a password depending on the security level set by the Supervisor

etc

Select any item to change




Use keys to select on or off




Press select to save any change or Menu to abort




Backlight

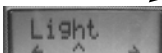


Vibrate Mode

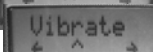


Beep on key press

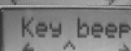

Sound on Alarm

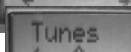

Sound if RF Fail

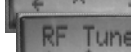

The RF Fail warning could be a normal daytime occurrence if the bedside units are turned off. Disabling the tunes stops the annoyance of the audible warning, but the display will still warn on the situation